

# STERIKING®



## **Rotosealer RS-120**

Economical choice for small clinics

- ▶ Simple operation
- ▶ All-day ready for use
- ▶ Increased sealing speed
- ▶ Constant sealing pressure
- ▶ Digital temperature control
- ▶ Sealing technique according to DIN 58953, PART 7
- ▶ Safety tested and approved by EN, VDE
- ▶ Built-in cooling control
- ▶ Practical worktop available as optional extra
- ▶ Shock watch system for delivery safety

**WIPAK MEDICAL**

## Rotosealer RS-120

continually ready for use for quick and reliable closing of all heat-sealable materials.

### Safe and quick sealing process


- Sealing temperature is easily adjustable to suit the material being sealed.
- When switched on, the motor and conveyor belts start automatically.
- The package is conveniently carried by the conveyor belts through the heating and sealing rollers to achieve quick and efficient closure.
- Adjustable in-feed guiding plate provides predefinition of the seal distance (0 - 35mm) from the package edge.
- Built-in fan improves the inside temperature control and minimizes overheating risk in all-day use.



The white epoxy coated housing is easy to maintain and to keep clean.



The package is hermetically sealed to protect against bacteria.

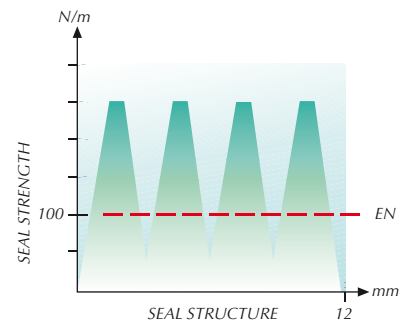

The Steriking® Rotosealers provide a strong seal. Total seal width is 12 mm.

Technical data	
Seal width	12 mm
Sealing speed	10 m/min
Sealing temperature	0 - 220 °C
Temp. tolerance	±2 %
Power consumption	250 W (max 600 W)
Power supply	220/240 V 50/60Hz
Weight	15 kg
Dimensions	490 x 255 x 160 mm
Separate worktop	625 x 225 x 60 mm

**Warranty**

- 24 months from the date of purchase.

**Shock watch label**

- Enables recognition of transport damage.
- Adhered to the transport carton.

Detailed operating manual and list of spare parts come with each unit.



**WIPAK MEDICAL**

Wipak Oy  
 P.O. Box 45, Wipaktie 2  
 FI-15561 Nastola, Finland  
 Tel. +358 (0)20 510 311  
 Fax +358 (0)20 510 3333  
 E-mail: steriking@wipak.com

www.steriking.fi www.wipak.com