


PercSys Accordion® Stone Management Device

Simplify
and Shorten
Endoscopic
Lithotripsy
Procedures




It's About Time!™

The PercSys Accordion Stone Management Device combines the utility of multiple endoscopic stone management tools into a single device, allowing stone localization and fragment removal using a single ureteroscopic instrument to simplify and shorten laser lithotripsy procedures.

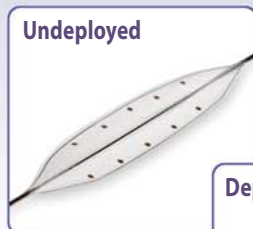


PercSys®
INNOVATIONS IN UROLOGY

PercSys Accordion® Stone Management Device

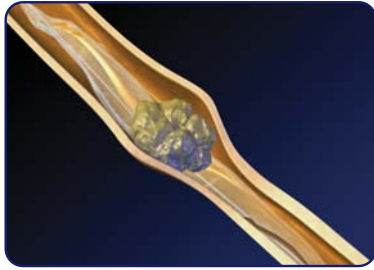


TRACKS like a guidewire

- Hydrophilic coating enables the device to track past impacted stones like a premium guide wire
- Flexible 6 cm tip facilitates entry and maintains access to the ureteral orifice
- Radiopaque markers allow easy monitoring of advancement and deployment using fluoroscopic or endoscopic visualization

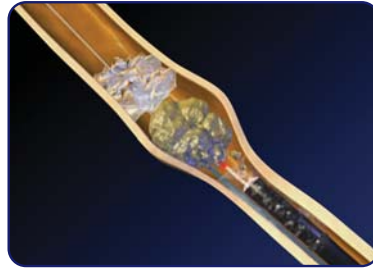


SECURES like a basket

- Multi-fold film occlusion is designed to prevent retrograde migration of the stone or small fragments
- Robust occlusion film withstands inadvertent contact with holmium laser and intracorporeal lithotripter energy
- Linear actuation enables the Accordion device occlusion to release easily from intact stones when required

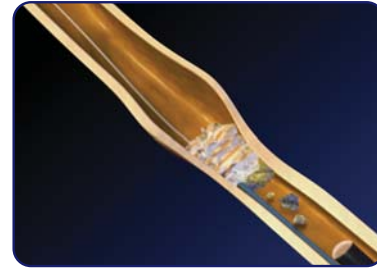


CLEARs like nothing else

- Film occlusion conforms to and fills the lumen of the ureter
- Allows increased irrigant inflow to improve visibility and continually rinse smaller lithotripsy fragments distally
- Sweeps larger fragments into the bladder or sheath following treatment

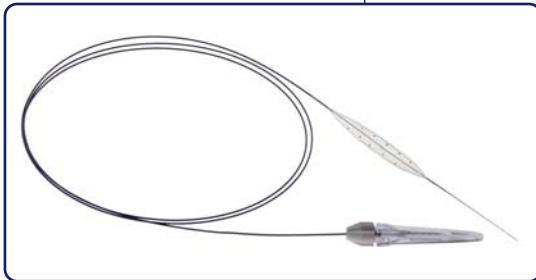

Undeployed

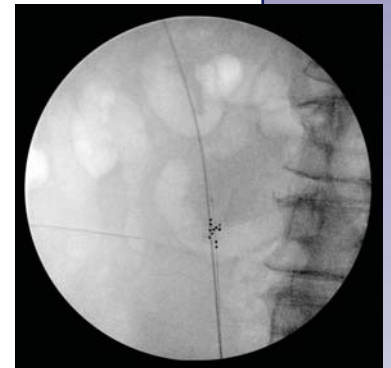

Deployed

Images courtesy of Demetrius Bagley M.D.
Thomas Jefferson University Hospital

Ordering Information

PercSys Accordion Stone Management Device

Deployed OD: (nominal)	Length:	Shaft OD:	Order No: †
7mm	150cm	0.038"	PA1205-06-07
10mm	150cm	0.038"	PA1205-06-10

† Available from Customer Service in boxes of 5 only.

Percutaneous Systems, Inc.
3260 Hillview Avenue, Suite 100
Palo Alto, CA 94304
Toll Free: 1.877.987.6633
Tel: 1.650.493.4200
Fax: 1.650.493.4212

www.percsys.com

PercSys
INNOVATIONS IN UROLOGY

© 2008 Percutaneous Systems, Inc. U.S. and foreign patents pending.
PercSys and PercSys Accordion are registered trademarks and *It's About Time!* is a trademark of Percutaneous Systems, Inc.
90-0750-01 Rev C 7/2008